

Biocraft's Pharmaceutical Equipments

- CPHM-8A **Craft's Collapsible tube Sealing Machine** Rs. 3,600
- CPHM-9 **Craft's Capsule Filling Machine**—Semi automatic machine for various operations with levers in order of (i) Putting of capsules with hands, (ii) Removal of Caps, (iii) Filling of Powder, (iv) Fixing of Caps and Sealing (v) Injection of capsules. It is for 100 capsules of size 0 to 2 or 5 Rs. 25,800
- CPHM-12 **Craft's Disintegration Test Machine**—Electrically operated, suitable for finding out disintegration tablets. Analog. Rs. 22,800
- CPHM-14 **Craft's Dissolution Rate Test**—For determination of dissolution rate test of tablet to USP STD complete with water constantly in motion and holding the temperature at 37 ± 0.5 deg. C through the test 1000 ml. Capacity cylindrical vessel with its slides are flanged near the top to accept a fitted cover having 4 parts. The motor speed is controlled with speed regulator in a range of 25 to 150 R.P.M. cylindrical stain less steel basket to hold the samples is of 316 quality single test. Rs. 19,500
- CPHM-15 **Craft's Tablet Friability Test Apparatus**—Special designed for finding out durability of tablets from the time of production to the time of use, as its important to test the abrasion and impact hardness of tablets—

—do—


single drum (Analog)

Rs. 7,900


—do—

but with single drum digital


Rs. 16,800

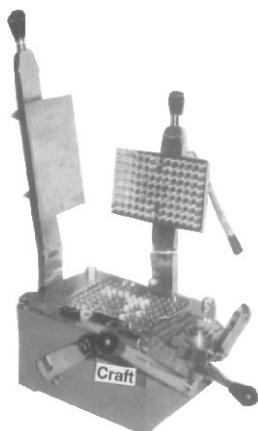


CPHM-9

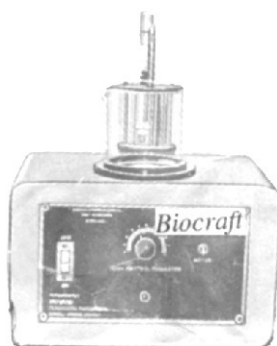

CPHM-12

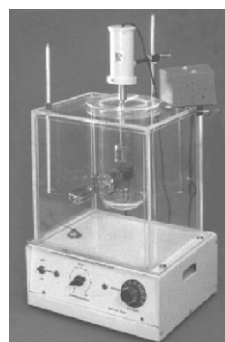

CPHM-14

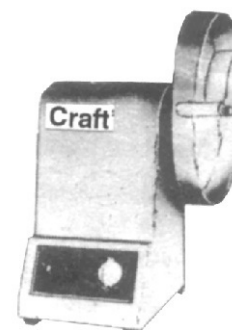


CPHM-15

- CPHM-16 **Craft's Filtration Equipment**—Having Following:
 (a) Membrane filter holder dia 47 mm. (b) Buckner Funnel of Glass,
 (c) Vacuum Pump and Flask. Rs. 14,800

- CPHM-19 **Craft's Hand Grind Mill**—With adjustment for making coarse and fine powder on metallic stand. Rs. 3,100

- CPHM-23 **Craft's Suppository Mould**—Made of brass for making
- | | | |
|----------------|----------------|------------------|
| 4 Moulds | 6 Moulds | 12 Moulds |
| 1,2,3,4 gm | 1,2,3,4 gm | 1,2,3,4,4 gm |
| Rs. 650 | Rs. 820 | Rs. 1,500 |


- CPHM-24 **Craft's Tablet Making Machine**—Size of the die and punch 9 mm. automatic compression and ejection device with revolution of the wheel drive. Automatic feed, removal of tablets. The unit is built on a heavy cast or mold steel fabricated with eccentric pressure system. Hopper provided and connected by a rubber tube. Easy setting for thickness and hardness of tablets. Complete with ONE Die and set of punch of chromium steel. Rs. 24,500

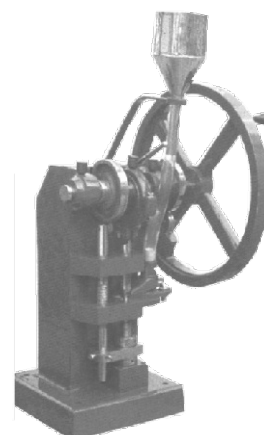

CPHM-24

Biocraft's Pharmaceutical Equipments


- CPHM-25 **Craft's Coating Pan 12" Made of Stainless Steel**—
 The machine fitted with gear box. Air blower is operated independently. An energy regulator is provided to control temperature. **Rs. 12,600**
- CPHM-26 **Craft's Tablets Polishing Pan**—Polishing pan made of mild steel provided with inner of superior canvas cloth 8" size **Rs. 3,600**
- CPHM-27 **Craft's Hardness Tester**—
 Monsanto type **Rs. 2,100**
 Pfizer type **Rs. 4,650**
- CPHM-28 **Craft's Tablet Counter**—For 100 tablets. **Rs. 1,100**
- CPHM-31 **Craft's Crown Sealing Machine**—Hand operated. **Rs. 6,100**
- CPHM-32 **Craft's Tincture press** **Rs. 3,600**
- CPHM-34 **Craft's Bottle Filling Machine-Motorized** single **Rs. 38,800**
- CPHM-35 **Craft's Bottle Filling Machine**—Electric operated MS body. **Rs. 11,800**
- CPHM-37 **Craft's Bottle Washing Machine** with 1/4 HP Motor—Single Brush **Rs. 9,500**
 Triple Brush **Rs. 13,700**
- CPHM-39 **Craft's Hot Extraction** (Soxhlet) equipment complete **Rs. 3,900**
- CPHM-41 **Craft's Bulk Density Apparatus** with Digital timer **Rs. 8,500**
- CPHM-45 **Craft's Single Stroke Motor Operated Tablet Machine**—Single stroke machine is an eccentric driven press for tablets upto maximum 1/2" diameter. The speed of the machine is 60 to 80 tablets per minute with single punch fitted. It can take up multi-punches as per the diameter of the tablet. Tablets which can be pressed may have round, square, triangular or any other design depending on the Die and Punch **Rs. 32,800**

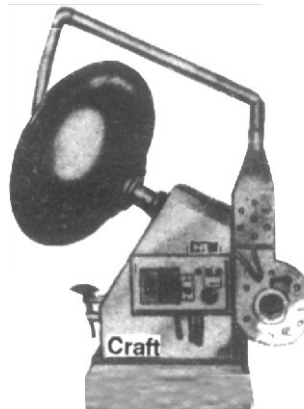


CPHM-25




CPHM-27

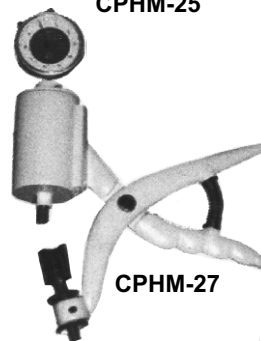

CPHM-37

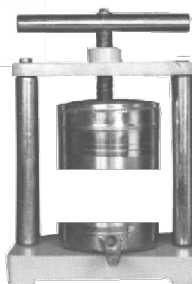

CPHM-32

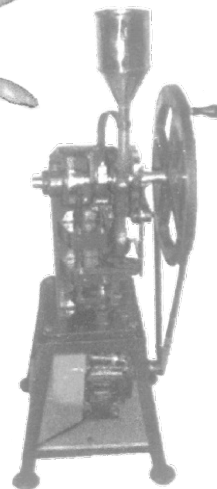

CPHM-45

Biocraft's Field Equipments

- BFE-2 **Field collection sack**—Made of heavy canvas with shoulder strap and pockets to hold 16 tubes—complete with tubes. **Rs. 570**
- BFE-3 **Prismatic Binocular**—Imported Routine quality **Rs. 1,600** each
- BFE-4 **Prismatic Binocular**—Imported double coated optics 7 × 35 or 8 × 40 or 10 × 50 (any) **Rs. 5,325** each
- BFE-5 **Prismatic Zoom Binocular** **Rs. 7,800** each



BFE-3



BFE-5